


iCluster User Manual v1.06beta

© 2016 Nick Hamilton

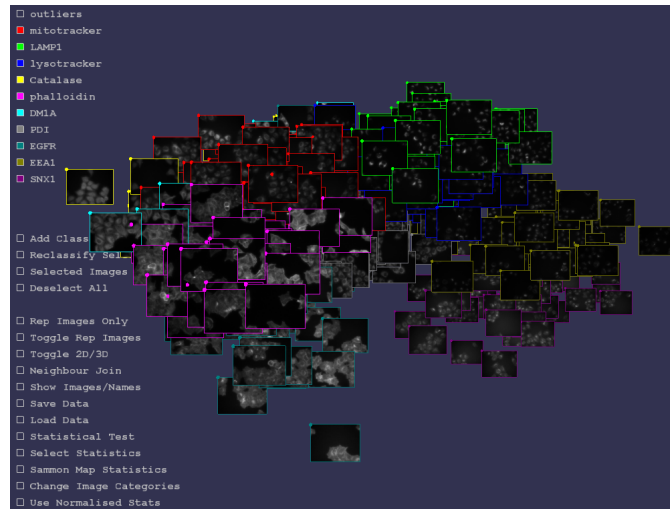

License

iCluster is released under Gnu General Public License. It comes with no warranty whatsoever. Details of the licence may be found in the licence file contained in the distribution.

Software © 2016 Nick Hamilton, n.hamilton@imb.uq.edu.au

Images © 2016 Rohan Teasdale, r.teasdale@imb.uq.edu.au

If you find this software useful in your research please reference paper (1) below.

Overview

The core idea of iCluster is to spatially place images of fluorescent protein subcellular localisation in such a way that images that are similar are spatially close. It does this by generating *threshold adjacency statistics* for each image to associate a 27 dimensional vector of real numbers with each image. The set of points in 27 dimensional space is then mapped into 2 or 3 dimensions in such a way as to preserve the distances between the points as well as is possible. Hence statistically similar images are spatially close in the visualisation. Statistical tests for difference between image sets may also be performed, statistically representative images found, along with a number of other features. The methods are described in [1,2,3]. Here, how to use the principle features of iCluster are described. Please email Nick Hamilton (n.hamilton@imb.uq.edu.au) with any comments, bug reports or requests.

Installation

Source code is in the subdirectory *src*. To compile the source code you will need the proce55ing software available from <http://proce55ing.org>.

Instructions for installation are contained in the README file in the *src* directory.

WARNING: iCluster has not yet been updated to work with Processing 3. Processing 2 is the preferred version to use.

Loading an Image Set

To load an image set, the images should be placed in a single directory. Images need to be 8-bit gray scale in either tiff or png format *if statistics are to be generated by iCluster* (any png image can be *displayed* in iCluster, but statistics generation requires 8-bit gray scale).

The “natural” format of iCluster is png, and tiff format images are automatically converted to png and saved to disk.

Images may be assigned to *classes*, for instance *treated* or *untreated*, in an image description file that is contained within the same directory as the image set. Classes are distinguished in the visualiser by border colours of the images visualised. The format for the file is given at the end of the document. An example image set together with an image description file (icluster_description.txt) can be found in the *images* folder in the distribution.

When iCluster starts, a file selection window will appear. Selecting an image description file will then load and display the image set.

Note: If images are missing, they will be displayed in iCluster as a black image with a red X.

Note: For large image sets, memory may become an issue. If iCluster detects that memory is running short, some images might not be displayed and a red X image will shown in their place. A warning message is also shown on screen which will suggest pressing the '-' key to automatically scale the images (and reduce memory usage) to a smaller size for display. It may be necessary to press '-' several times until the warning message goes away.

Note: If an image set is loaded for the first time, statistics will be automatically generated for that set if statistics have not been supplied in the image description file. However, the initial placement of the images in the viewing environment will be random. To find an optimal spatial placement of the images so that those that are statistically similar are spatially close click on the Sammon Map check box (see below).

Navigating an Image Set

The following mouse and keyboard controls may be used to navigate around an image set:

View/Hide an Image Class	Left mouse click on the check box for that class
Zoom In/Out	Mouse scroll
Rotate/Translate View	Left mouse drag
Increase/Decrease Image Size	Right mouse + mouse scroll + CTRL
Increase/Decrease Sphere Size	Right mouse + mouse scroll + SHIFT
View image in detail	Left mouse on image
Select an Image	Right mouse on image
Select a Group of Images	Right Mouse on a central image and drag to select around it
Randomise Image Positions	/
Start/Stop Sammon Mapping	.
Show/Hide Image Masks	m
Increase/Decrease Image Resolution	+/-
Snapshot Image (to Sketch Folder)	p

Note: for memory reasons, images are initially viewed at 1/3 of their size. This resolution may be increased/decreased using +/- . Similarly, displaying image masks (see below) doubles the amount of ram used and may cause problems for large image sets at high resolution. A recent dedicated graphic card is recommended for use with iCluster to get the best performance.

2D/3D Representation

Image sets may be viewed layed out in 2D or 3D. 3D will typically be better in the sense that the Sammon mapping (see below) of statistics will better preserve the distance relationships between statistics for the images. Clicking the Toggle 2D/3D button switches between these representations.

Show Images/Names

The data may either be viewed as a set of images spatially placed by similarity, a set of names of those images similarly placed, or as a sphere. Clicking on 'Show Images/Names' toggles between these.

Nearest Neighbours

The nearest neighbour to an image is the image that is 'nearest' to that image. Nearest may either be in the sense of the 2D or 3D distance in the viewer, or in the (high dimensional) Euclidean distance between statistics vectors. Clicking in the 'Neighbour Join' will toggle between showing lines joining nearest neighbours in each of these senses. The number of images whose nearest neighbours are the same class is also displayed. This gives a sense of how separated the images classes are.

Note: Results may vary according to whether 'Normalise Statistics' (below) and 'Select Statistics' (above) are set.

From iCluster v1.04 *consecutive* images may also be joined. That is, images that are consecutively described in the iCluster description file are joined (if they are of the same class). The idea here is that the frames of one or more time lapse movies can be visualised with the consecutive frames visually connected by lines.

Representative Images

A representative image for each class may be shown either using the 'Rep Images Only' checkbox (where images are shown side by side) or 'Toggle Rep Images' where representative images are shown spatially located in 2D or 3D and all other images are hidden. A representative image is chosen for each class as that image which is nearest to the centroid of the image statistics for that class. Image size may be increased/decrease in this mode using Right Mouse Button + Scroll Mouse Button.

Note: Results may vary according to whether 'Normalise Statistics' (below) and 'Select Statistics' (below) are set.

Add Class

New classes to which images may belong may be added by clicking on the 'Add Class' checkbox. When selected, the user is prompted to enter a name for the class. Press 'Enter' to conclude entering the name. One use for this is to create an 'Outlier' class of images which are unusual in some way.

Reclassify Selected

Images that have been selected via right mouse clicking appear with a red tint. Such a set may then be reclassified to another class by clicking the 'Reclassify Selected' checkbox. Once checked, then user can choose which class to reclassify the images to by clicking on one of the class checkboxes.

Deselect All

All of the selected images may be deselected by clicking the 'Deselect All' checkbox.

Selected Images Only

If images have been selected via right mouse clicking, that set of images may be viewed side by side in 2D by clicking the 'Selected Images Only' checkbox. Exit this mode by clicking the

checkbox again. Image size may be increased/decrease in this mode using Right Mouse Button + Scroll Mouse Button.

Load/Save Data

An image set may be loaded into iCluster by clicking the 'Load Data' button. By selecting either a tiff or png format image in a directory using the file browser that appears, all of the images in that directory will be loaded into iCluster. However, images so loaded will all be of the 'Default' class.

To load images with classes also supplied, a Image Description File (described below) must be in the same directory as the images. Selecting this file will load the images, their associated classes, and any statistics or coordinates previously calculated.

An image description file may be saved by clicking on 'Save Data' and selecting or creating a file *in the same directory as the images*. The next time iCluster is started the last saved Image Description File will be automatically loaded along with the associated images.

Statistical Testing

Classes of images may be tested against each other to detect if they are statistically different. When the 'Statistical Test' checkbox is clicked, the user is presented with 2 lists of the image classes. Those image classes selected in the left hand list are then tested against those selected in the right hand list. Once classes have been selected, click the 'Finished Selecting Classes' checkbox, and the test will start.

The test is for the null hypothesis, that the image sets come from the same population, more specifically that the mean vectors of the population that the two sets are sample from are the same. This is performed using permutation testing as described in [1]. A p-value for the null hypothesis is returned. The number of repeats to give the p-value may also be selected, though under most circumstances, the default of 1000 should be sufficient. The Euclidean distance between the centroids (average statistics vector) of each cluster are also returned.

To exit the Statistical Testing mode, click the 'Statistical Test' checkbox.

Note: Results may vary according to whether 'Normalise Statistics' (below) and 'Select Statistics' (below) are set.

Selecting Statistics

The statistics that are used within iCluster may be selected or deselected by clicking on the 'Select Statistics' checkbox, and then checking/unchecking the listed statistics. It is probably safe to ignore this unless you want to use iCluster with your own image statistics rather than the automatically generated *threshold adjacency statistics*. If you do use your own image statistics within iCluster, the ability to select/de-select subsets of them may be useful if you wish to examine the effect on clustering or statistical testing of just using certain of your statistics. See 'Sample Image Description File' below for examples of how to use your own statistics in iCluster.

Image Categories

This is probably not of interest to most users, but iCluster as well as supporting having a class associated with each image, also supports having categories of classes. For instance the category 'Subcellular localisation' might contain the classes 'Nucleus', 'Cytoplasm', 'Mitochondria' and so on. Most image sets will have just one category containing the classes. However, it is possible to have other categories of classes, and you can switch between these by clicking on the 'Change Image Categories' checkbox. Another category might be one that contains the classes 'Punctate' and 'Not Punctate'. Hence it is possible to have multiple categories with a class from each category associated with each image, and then to switch between views of each category. Categories and the classes within those categories may be defined in the Image Description File, and an example may be found at the end of this document.

Sammon Mapping

Sammon mapping is the process by which high dimensional statistics vectors are mapped into 2 or 3 dimensions in such a way as to preserve distances between points as well as is possible. Click on 'Sammon Map Statistics' to start the process. It is iterative, and so you will see the images fly around the screen a bit as they find their preferred location. After a period of time the images should become relatively stable in their location, at which point you may click the checkbox again to stop the process.

Note: Results may vary according to whether 'Normalise Statistics' (below) and 'Select Statistics' (above) are set.

Principal Component Analysis (PCA)

Images may be positioned in 2D or 3D using the principal components of the statistics vectors for the images. This may be useful to observe trends in the data, but more importantly the spatial layout using PCA can be used as an initial configuration for Sammon mapping (rather than starting with randomly positioned points). In this way results of Sammon mapping will be reproducible. Also, testing has shown that the Sammon mapping converges more quickly to a stable configuration when initialised with PCA coordinates.

Normalising Statistics

Statistics may be normalised, that is for each statistic then mean of the statistic is subtracted and the result divided by the standard deviation for that statistic. Threshold adjacency statistics are *standardised* for the complete image set when first calculated (see below). However, if you wish to only compare some subset of the complete image set it may be appropriate to normalise the statistics with regard to that subset. Hence when "Normalise Statistics" is checked, normalisation is performed with regard to those classes of image that are *visible*, i.e. their class checkboxes are selected. Calculations such as representative image selection, high dimensional nearest neighbours, sammon mapping are then performed with respect to the normalised statistics for the visible images. In the case of statistical testing, normalisation occurs with regard to those classes of images that have been selected for comparison.

Threshold adjacency statistics are *standardised* when first calculated as follows. The normalisation means and standard deviations of the TAS statistics for the 500 images of 10 subcellular localisations in Image Set A (described in [1]) were calculated. For any new image set, the TAS as calculated, and then each statistic is *standardised* by subtracting the mean and dividing by the standard deviation for that statistic *as calculated for the images for Image Set A*. The idea is that Image Set A provides a reference by which to standardise the statistics for any image set. Hence when standardised, the Euclidean distances between images statistics vectors are comparable across any image sets for which the statistics have been standardised. In other words, standardisation means we have a metric to measure the distance between any images. To use normalised rather than standardised statistics select the "Normalise Statistics" check box.

Show/Hide Image Masks ('m' key)

Statistics are generated for each image by first selecting the cellular regions in the image using a simple thresholding scheme. The threshold is calculated as $m + 0.9\sigma$, where m is the modal intensity in the image, that is the most common intensity value, and σ is the standard deviation of the intensities in the image. While this produces a reasonable selection for the majority of image iCluster has been tested on, it may not work in all cases such as when the background is uneven. In such a case, the result is often an outlier image.

Masks may be viewed by pressing the 'm' key while images are being viewed. Press 'm' again to

return to viewing the images.

In cases where the automated masking fails completely, the user may supply their own masks for each image. A mask for an image is an image for which the non-zero intensity pixels select the regions of interest in the original image. See Sample Image Description File section below for how to supply a mask image for an image.

Snapshot Images

Pressing 'p' will take a snapshot image of the iCluster window. The image is automatically saved to the directory containing the iCluster executable.

Sample Image Description Files

There are a number of formats which will allow images/classes/statistics to be defined in an iCluster Image Description File. Essentially the file is plain text with a single line for each image, comma separated, describing its name, class(es), possibly a mask, possibly user supplied statistics. Hence exporting an Excel document with appropriate column entries to comma separated format text file (csv) should produce a correctly formatted file, though the column headings may need to be removed.

In it's simplest form, in each line an image name and a class is given, for instance:

```
010_mitotracker.png,mitochondria
```

defines the image 010_mitotracker.png to belong to the class 'mitochondria'.

To supply a mask with each image, a each line might look like:

```
010_mitotracker.png, 010_mitotracker_mask.png,mitochondria
```

which is as before, but with '010_mitotracker_mask.png' as the second entry, that being the name of the mask image associated with the image 010_mitotracker.png. The fact that the second entry is a mask image is recognized by it's name terminates in either '.png' or '.tif'.

To supply your own statistics for each image, a line might look like

```
010_mitotracker.png,mitochondria,0.12,0.2,1.45,11.2,2,1.7,-0.3,7.2,-1.4
```

where the sequence of comma separated numbers are the statistics to use.

Note: Missing values are not supported by iCluster, so ensure all lines contain the same number of comma separated real numbers.

To have multiple categories for the classes of an image (see Image Categories, above) a list of classes may be given, instead of a single class. For instance:

```
010_mitotracker.png,mitochondria,not punctate
```

might describe the image 010_mitotracker.png as being of mitochondria (in the subcellular localisation category), but also being of class 'not punctate' (in the category punctuate/not punctate). Similarly, user supplied image statistics might be appended to such an entry.

The most important point for the Image Description File is that every line *must have exactly the same number of comma separated entries*. Hence if a mask image is supplied for one image, masks must be supplied for all. Similarly, each image must have the same number of image statistics (if supplied), and if there are multiple categories, each containing an image class, then all images must have entries defining what their class is within that category.

References

1. Hamilton N., Wang J., Kerr M.C., Teasdale R.D. Statistical and visual differentiation of high throughput subcellular imaging. BMC Bioinformatics 2009;10:94
2. Hamilton N., Pantelic R., Hanson K., Karunaratne S., Teasdale R.D. Fast automated cell phenotype image classification. BMC Bioinformatics 2007;8:113.
3. Hamilton N., Teasdale R.D. Visualizing and clustering high throughput sub-cellular localization imaging BMC Bioinformatics 2008;9:81.

Reference (1) describes the main methods and testing of iCluster. This is the article to reference if you find iCluster useful in your research.

Reference (2) describes the *threshold adjacency statistics* that are utilised by iCluster for visualisation and statistical testing.

Reference (3) describes the method of mapping statistics into 2 or 3 dimensions for visualisation and an early prototype of iCluster.

iCluster is built using the fantastic *processing* programming language.

See <http://www.processing.org> .

Problems

Text entry on some Mac notebooks

A bug has been observed on a couple of Mac notebooks whereby text entry does not work, for instance, when adding a new class name. The source of the bug has not been tracked down yet, but there is a simple work around: after starting iCluster, minimise then then un-minimise the iCluster window.

iCluster runs but is sluggish when viewing image

A limiting factor in the visualisation is the graphics card of the computer. A poor graphics card, or integrated graphics, may produce poor results for larger image sets. A better graphics card will improve things. Alternatively, the resolution that images are held in memory at may be reduced by pressing '-' while iCluster is running. This will reduce the amount of memory needed for the visualisation, but will mean that images are not seen in such detail. iCluster has successfully viewed a 1000 images of individual cells simultaneously by using a good graphics card on a notebook computer and reducing the image resolution.

iCluster doesn't do anything when starting

Try deleting the file config.txt in the directory with the executable. iCluster may be getting confused about where the default directory/image set is.

Again, large image sets may cause problems due to memory limitations. Windows Vista seems to be particularly inefficient in this respect. The Linux distribution of iCluster will readily display twice the number of images that Vista will on identical hardware.

One work around for this is to reduce the image resolution displayed when iCluster starts. In the folder containing the executable is the file *config.txt*. This file contains 3 lines of information, the default directory of images, the name of the image description file, and a line that reads

iScale=n

where n is an integer. This is a scaling factor used to reduce the image size. For instance, if n=3, the images will be scaled to 1/3 of their original size. Hence larger image sets may be displayed by increasing n. Note that this only effect the display of the images, statistics are still calculated on the raw images.

The input parsing of the description file currently does little checking that the input is well formed. Hence if you have created your own description file, one potential source of iCluster freezing is a poorly formed input file with missing entries and so on. Later versions of iCluster will check that input is well formed and issue appropriate error messages.

Additional Coding:

PCA implemented by Daniel Marshall using the jmathtools library

<http://jmathtools.berlios.de/doku.php>